

dosasap®

PC controlled dosimetry equipment for CTA dosimeter readout

# A COMPACT TOOL FOR IQ/OQ/PQ PURPOSES




## SOFTWARE PACKAGE MAIN FEATURRES

- Dosimetry system calibration (readout of calibration CTA strips, curve fitting, residual and uncertainty calculations)
- Electron Beam width, scan length and scanning uniformity assessment
- Photon and Electron Beam absorbed dose profile measurement for dose uniformity estimation in user selected Region Of Interests (ROI) strip length from 5 cm to 12 m
- Electron Beam energy determination with wedge technique according to ISO/ASTM 51649

aerasap®

Software for continuous optical dosimetry tools

data base driven compliant with relevant standards ISO/ASTM 51649, 51650, 51261, ISO 11137-3, ... and FDA CFR 21 part 11.

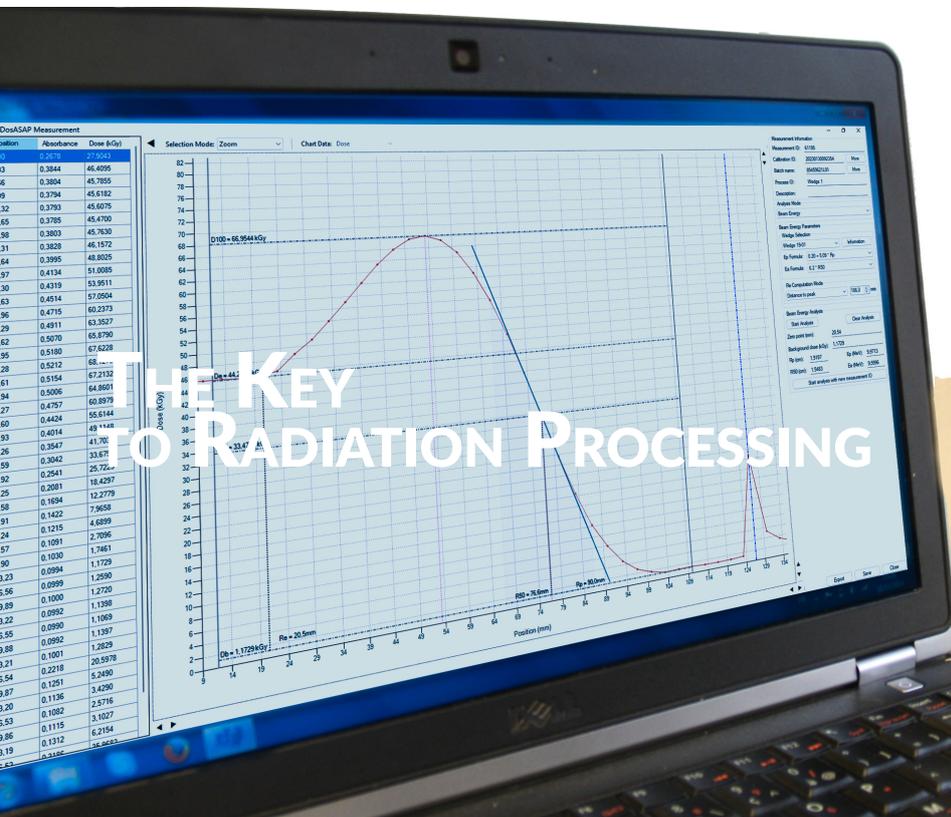

# THE KEY TO RADIATION PROCESSING

## Contact us

- Sébastien Riegler - Product Manager - Sales - After Sales  
s.riegler@aerial-crt.com
  - Abbas Nasreddine - Project Manager - Development  
a.nasreddine@aerial-crt.com
- +33 (0)3 88 19 15 15  
Parc d'innovation - 250 rue Laurent Fries 67 400 Illkirch - France  
[www.aerial-crt.com](http://www.aerial-crt.com)

